


1 kb DNA Ladder

Ready to load agarose size standard

Containing 14 double stranded DNA fragments (34 ng/μl) in a room temperature stable loading buffer

Description

Fragment sizes from top to bottom = 10000, 8000, 6000, 5000, 4000, 3500, **3000**, 2500, **2000**, 1500, 1000, 750, 500 & 250 bp

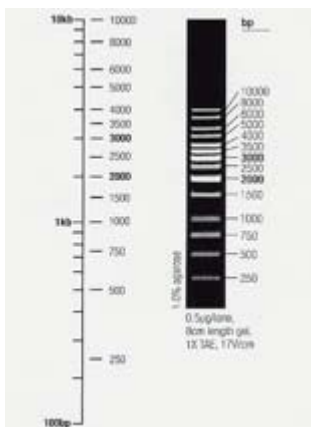

A genuine Microzone product, made in the UK

Load 10 to 15 μl per well in the agarose gel, alongside samples

Related Products

MegaMixes = Clear or blue PCR mixes for instant and accurate PCRs

100 bp DNA ladder = An agarose size standard supplied 'ready to load' in room temp. stable loading buffer

microCLEAN = The easy way to clean-up PCR fragments

Supplied as a blue solution in multiple of 1 ml aliquots
Store at +4°C or at room temperature

For Research Only